


Guardian In-line Particle Filter

Eliminate blockages

If there are particulates in your samples, there is a risk that they may get trapped in the fine channels of your sample line or within the nebulizer. The Inline Filter provides a simple and effective way to eliminate this risk. It incorporates a 120 micron filter and is easily inserted in the sample tubing between the autosampler probe and the nebulizer.

The purpose-built clog-resistant design is ideal for ICP samples. Any particle build-up is easily removed by back-flushing and ultrasonic cleaning. And the PEEK material is suitable for use with all of the most common ICP solutions.

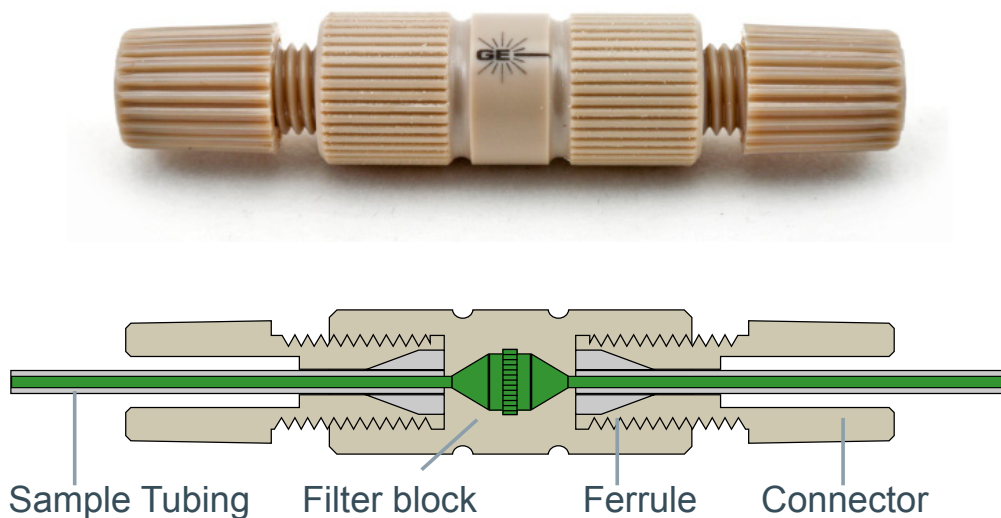

For more information visit: www.geicp.com or
contact us at enquiries@geicp.com



GLASS EXPANSION
Quality By Design

Asia Pacific

6 Central Boulevard
Port Melbourne VIC 3207
Australia
+61 3 9320 1111
enquiries@geicp.com

Americas

31 Jonathan Bourne Drive
Unit 7, Pocasset, MA 02559
USA
508 563 1800
geusa@geicp.com

Europe

Friedenbachstrasse 9
35781 Weilburg
Germany
+49 6471 3778517
gegmbh@geicp.com